


NP 15

The Quietest Nasal Pillow Mask for High-Level Safety and Wearing Comfort

Only 15 dB at 10 hPa – our nasal pillow mask NP 15 offers high performance in all areas. It is quiet, small and yet stable and efficient. Moreover, the easy-to-use NP 15 features an almost unlimited range of anchoring options to keep the tube securely in place. Patients enjoy greater wearing comfort and consequently improve their therapy compliance. The ball-and-socket joints, one on the mask body and the other between the mask and tube, allow the patient unrestricted freedom of movement while ensuring a secure fit.

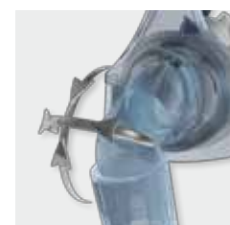

The two ball-and-socket joints allow...



...unrestricted freedom of movement

NP 15, WM 26755

Therapy Comfort

- The quietest nasal pillow mask
- Has two ball-and-socket joints:
 - The ball-and-socket joint on the mask body allows unrestricted freedom of movement.
 - The second ball-and-socket joint lies between the mask and device tube and ensures a stable fit free of leaks.
- Has two anchoring options (clip and strap) which help to hold the tube reliably in place throughout the night.
- Integrated pressure measurement port
- No annoying drafts. The exhaled flow of air is directed away from the face and the outflow gap toward the face is closed.
- Optimum fit
- Small and delicate but stable and efficient
- Made of skin-friendly silicone

Operating Convenience

- The different colors of the Velcro tabs (light blue above, dark blue below) make the headgear easy to use.
- Pillows move independent of each other
- The crossband on the headgear is adjustable.
- Three pillow sizes
- Vertical adjustment on forehead support makes mask fitting fast and simple
- Immediately ready for use; all essential parts are in standard delivery package
- Adjustment range clearly marked on forehead support
- Small number of parts keeps operation simple
- Headgear lies above the ears, making it easier to put on and take off the mask.
- The many replacement parts help to reduce costs.

Device Description	WM No.
NP 15	26755



Supply Schedule

- 1 Nasal pillows, Size S (not shown) – WM 26782
Nasal pillows, Size M (not shown) – WM 26787
Nasal pillows, Size L – WM 26792

- 2 Headgear including
Tube anchoring options – WM 26775

- 3 Forehead support set – WM 15623

- 4 Forehead cushion – WM 26777

- 5 Mask tube set – WM 15744

- 6 Tube anchoring clip – WM 26728

Tube anchoring strap (not shown) – WM 26774

Set of headgear clips (2 pieces) (not shown)
– WM 15622

Our complete range of therapy solutions, accessories and mask systems is at: weinmann.de

Technical Data NP 15



Product class as per 93/42/EWG:	II a	Temperature range	
Dimensions (W x H x D):	approx. 140 x 100 x 40 mm	■ Operation:	+5 °C to +40 °C
Weight:	about 90 g (depending on size and configuration)	■ Storage:	-20 °C to +70 °C
Therapy pressure range:	4 hPa to 25 hPa	Flow resistance	
Tube connection:	Cone as per DIN EN ISO 5356-1 Ø 22 mm (male)	■ at 50 l/m:	≤ 0.9 hPa
Pressure measurement port:	Ø 4 mm	■ at 100 l/m:	≤ 3.5 hPa
		Sound level	
		■ at 10 hPa:	15 dB(A)
		Service life:	up to 12 months
		(Materials used in the manufacture of the masks age and wear out when treated with aggressive cleansers. In isolated cases, therefore, it may be necessary to replace masks more frequently.)	

83548-EN-02/2013 © Copyright Weinmann, Hamburg. Duplication of any kind only with the express permission of Weinmann. All rights to design and specification modifications reserved.